


**INTEGRATED BLOOD CIRCUIT FOR SIMPLICITY
– IN JUST ONE MOVE**

As a pioneer and leading innovator in dialysis therapy, we present the **Artiset** blood tubing system that raise treatment quality, staff satisfaction and operational efficiency compared to conventional blood tubing systems.

The **Artiset** blood tubing system is an exclusive blood circuit designed to work with the **Artis Physio Plus** system. Compared with conventional bloodlines, the compact and highly integrated design offers various important advantages both from a medical and handling perspective.

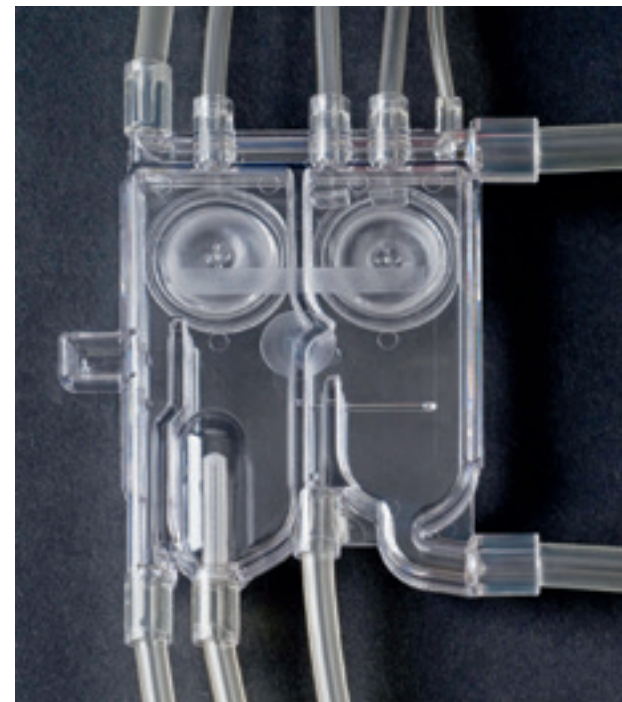

Baxter
Artiset
BLOODLINE TUBING SYSTEM



HD DN – ARTISET BLOOD TUBING SYSTEM

TREATMENT MODALITY	CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
Hemodialysis – Double Needle Saline Priming	955075 955549	Artiset HD DNL HC	Gamma/Beta radiation Beta radiation	20
Hemodialysis – Double Needle Ultra Prime	955075 955549	Artiset HD DNL HC	Gamma/Beta radiation Beta radiation	20
Hemodialysis – Double Needle Ultra Prime	955077	Artiset Prepost	Gamma radiation	20

ACCESSORY

CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
n.a.	n.a.	n.a.	n.a.
955600	Ultraline HD	Beta radiation	80
955600	Ultraline HD	Beta radiation	80

HD SN – ARTISET BLOOD TUBING SYSTEM

TREATMENT MODALITY	CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
Hemodialysis – Single Needle Saline Priming	114533	Artiset HD SN HC	Gamma radiation	7
Hemodialysis – Single Needle Ultra Prime	114533	Artiset HD SN HC	Gamma radiation	7

ACCESSORY

CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
n.a.	n.a.	n.a.	n.a.
955600	Ultraline HD	Beta radiation	80

ON-LINE HDF/HF IN PRE/POST DILUTION – ARTISET BLOOD TUBING SYSTEM

TREATMENT MODALITY	CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
Hemodiafiltration pre and postdilution	955077	Artiset Prepost	Gamma radiation	20
Hemofiltration predilution	955077	Artiset Prepost	Gamma radiation	20
Hemodiafiltration postdilution	955075 955549	Artiset HD DNL HC	Gamma/Beta radiation Beta radiation	20

ACCESSORY

CODE	DESCRIPTION	STERILIZATION	PIECE/BOX
955599	Ultraline HDF	Beta radiation	40
955599	Ultraline HDF	Beta radiation	40
955599	Ultraline HDF	Beta radiation	40

Note: A contingency stock of saline priming bags is recommended in case **Ultra** Prime function cannot be used.

RECOMMENDED SPARE ACCESSORIES

INTENDED USE	CODE	ACCESSORY	STERILIZATION	PIECE/BOX
Connector to pause treatment	114376	Evaclean rinse-backconnector	Ethylene Oxide	150
Infusion line for saline priming with Artiset Prepost	107472	Infusion line SP-339G	Beta radiation	60

For the safe and proper use of the devices referenced within, refer to the complete Instructions for Use.

CE 0086

Baxter, Artis Physio Plus, Artiset, Hemocontrol and Hemoscan are trademarks of Baxter International Inc. or its subsidiaries.

GBU-RC7-200002 v1.1 – October 2020

renalcare.baxter.com
Baxter Healthcare Corporation
One Baxter Parkway
Deerfield, IL 60015
USA
1-800-422-9837


ARTISET BLOODLINES
for the **Artis Physio Plus** system

HEMODIALYSIS

Hemodialysis Double Needle Artiset HD DNL HC

The Artiset HD DNL HC blood tubing system for double-needle hemodialysis treatment includes:


- The **Hemoscan** cuvette to measure blood volume variations, enabling the **HemoControl** modality. For double-needle hemodialysis only
- A connector on the venous dialyzer line for **Ultra Prime** function
- A pre-connected infusion line for saline priming

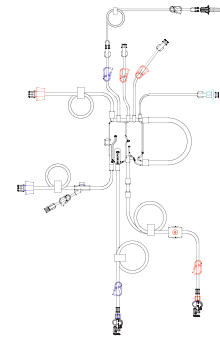

Hemodialysis Single Needle Artiset HD SN HC

The Artiset HD SN HC blood tubing system for single-needle hemodialysis treatment includes:

- The **Hemoscan** cuvette to measure blood volume variations
- The arterio-venous expansion chambers for effective single-needle mode
- A connector on the venous dialyzer line for **Ultra Prime** function
- A pre-connected infusion line for saline priming

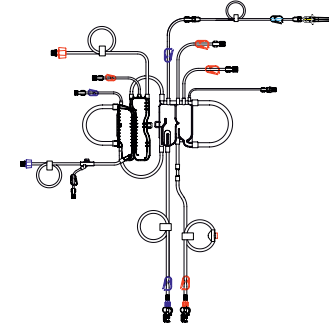


On-line auto-priming Ultraline HD

The one-button priming can be performed with on-line produced fluid using the **Ultraline HD** simply connected to the Ultra port. The **Ultraline HD** is intended for priming, bolus and rinseback in Hemodialysis (HD) Double Needle (DN) and Single Needle (SN) treatments with infusion fluid delivered on-line by **Artis Physio Plus** dialysis machine. Storage and transportation of priming bags is no more necessary.

This bloodline accessory is dedicated for the **Ultra Prime** function in hemodialysis treatments, to be used in combination with the **Artiset Bloodline** range and include:

- The T-connector for rinse-back to facilitate the rinse-back operation

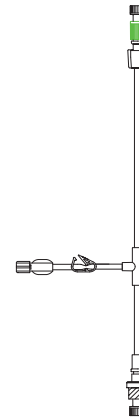

ARTISET BLOOD TUBING SYSTEM FOR HD TREATMENT

PRODUCT DESCRIPTION	ARTISET HD DNL HC	ARTISET HD SN HC
Product code:	955075, 955549	114533
Volume [ml]:	132	227
Shelf life:	3 years	3 years
Sterilization:	Gamma/Beta radiation	Gamma radiation
Units per box:	20	7

MATERIALS

Tubing:	PVC polyvinyl chloride (medical grade, DEHP-free)
Cassette body:	PETG polyEthylene Terephthalate Glycol
Injection site:	Latex-free rubber

ACCESSORY

PRODUCT DESCRIPTION	ULTRALINE HD
Product code:	955600
Volume [ml]:	n.a.
Shelf life:	3 years
Sterilization:	Beta radiation
Units per box:	80

MATERIALS

Tubing:	PVC polyvinyl chloride (medical grade, DEHP-free)
---------	---

PERFORMANCES


Max. infusion vol.:	15 l.
---------------------	-------

ON-LINE HDF – HF IN PRE/POST DILUTION

HDF pre or postdilution HF predilution Artiset Prepost

The Artiset Prepost blood tubing system for double-needle hemodialysis, HDF pre and postdilution, or HF predilution treatment includes:


- The **Hemoscan** cuvette to measure blood volume variations, enabling the **HemoControl** modality
- A connector on the venous dialyzer line for **Ultra Prime** function

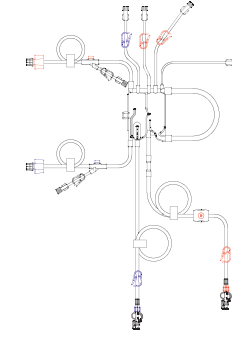

HDF postdilution Artiset HD DNL HC

The Artiset HD DNL HC blood tubing system for double-needle hemodialysis, HDF post dilution treatment includes:

- The **Hemoscan** cuvette to measure blood volume variations, enabling the **HemoControl** modality. For double-needle hemodialysis only
- A connector on the venous dialyzer line for **Ultra Prime** function
- A pre-connected infusion line for saline priming

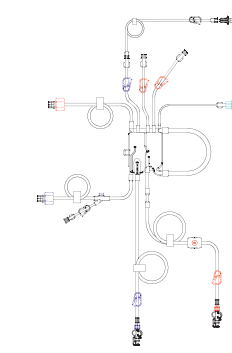


ON-LINE substitution fluid Ultraline HDF

The **Ultraline HDF** for Haemodiafiltration (HDF) Double Needle treatment in pre or post dilution modality and Hemofiltration (HF) Double Needle treatment in pre dilution modality together with on-line priming, bolus and rinseback of the blood circuit using infusion/substitution fluid delivered by **Artis Physio Plus** dialysis machine.

This bloodline accessory is dedicated for on-line modalities, to be used in combination with the **Artiset** range and includes:

- A combi-connector for auto-loading and unloading operation
- The T-connector for rinse-back to facilitate the rinse-back operation

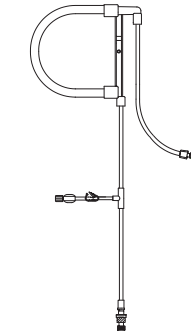

ARTISET BLOOD TUBING SYSTEM FOR HDF AND HF TREATMENT

PRODUCT DESCRIPTION	ARTISET HD DNL HC/ARTISET PREPOST
Product code:	955075, 955549/955077
Volume [ml]:	132
Shelf life:	3 years
Sterilization:	Gamma or Beta Gamma radiation
Units per box:	20

MATERIALS

Tubing:	PVC polyvinyl chloride (medical grade, DEHP-free)
Cassette body:	PETG polyEthylene Terephthalate Glycol
Injection site:	Latex-free rubber

ACCESSORY

PRODUCT DESCRIPTION	ULTRALINE HDF
Product code:	955599
Volume [ml]:	n.a.
Shelf life:	3 years
Sterilization:	Beta radiation
Units per box:	40

MATERIALS

Tubing:	PVC polyvinyl chloride (medical grade, DEHP-free)
---------	---

PERFORMANCES

Max. infusion vol.:	150 l.
---------------------	--------

